

OV5 HOMOGENIZER



WHAT IS DISPERSING?

The term “dispersing” means the **crushing** and **dispersing** of a solid, liquid or gaseous phase in a fluid with **which it cannot be totally mixed**.

- Dispersing is used to reduce the particle size of a material (crushing)
- Dispersing is also used to mix materials that are usually immiscible



DISPERSING

TERMS USED IN DISPERSING TECHNOLOGY

SUSPENSION - A coarse dispersing system of insoluble particles of solids in a liquid dispersion medium.

Example: X-ray contrast media, clay slurries, milk of lime, pharmaceutical preparations.

EMULSION - A thermodynamically unstable system consisting of at least 2 liquids which are normally immiscible, of which one is dispersed in spherules in the other liquid phase.

Example: lotions, creams, milk, soluble oil, margarine, mayonnaise.

AEROSOL - a dispersion of gas in liquid.

Example: Shaving foam



VISCOSITY

One of the most important criteria when working with dispersers is the viscosity. Viscosity is a measure for the internal friction or flow resistance of a liquid.

As viscosity increases, the flow rate decreases. Generally, **for each 1000 mPas increase in viscosity, the flow-rate decreases by 20-25%.**

If the viscosity or flow limit is too high, the material does not flow, which might lead to a flow interruption in the rotor-stator system.

Some examples of viscosity for common substances at 20 °C, in mPas:

- Water: 1
- Ethyl Alcohol: 1.2
- 100% Glycerine: 1,500
- Honey: 10,000


VELP OV5 is able to process viscosities up to 10'000 mPas



OV5 – ROTOR-STATOR PRINCIPLE

HOW IS FLUID SHEAR INDUCED IN ROTOR/ STATOR?

- Rotor rotates on shaft axis enclosed in stator
- Fluid enters on the same axis as drive shaft
- Momentum is changed 90 degrees and moves in radial direction due to centrifugal force
- Fluid enters channel between rotor teeth subjected to centrifugal thrust

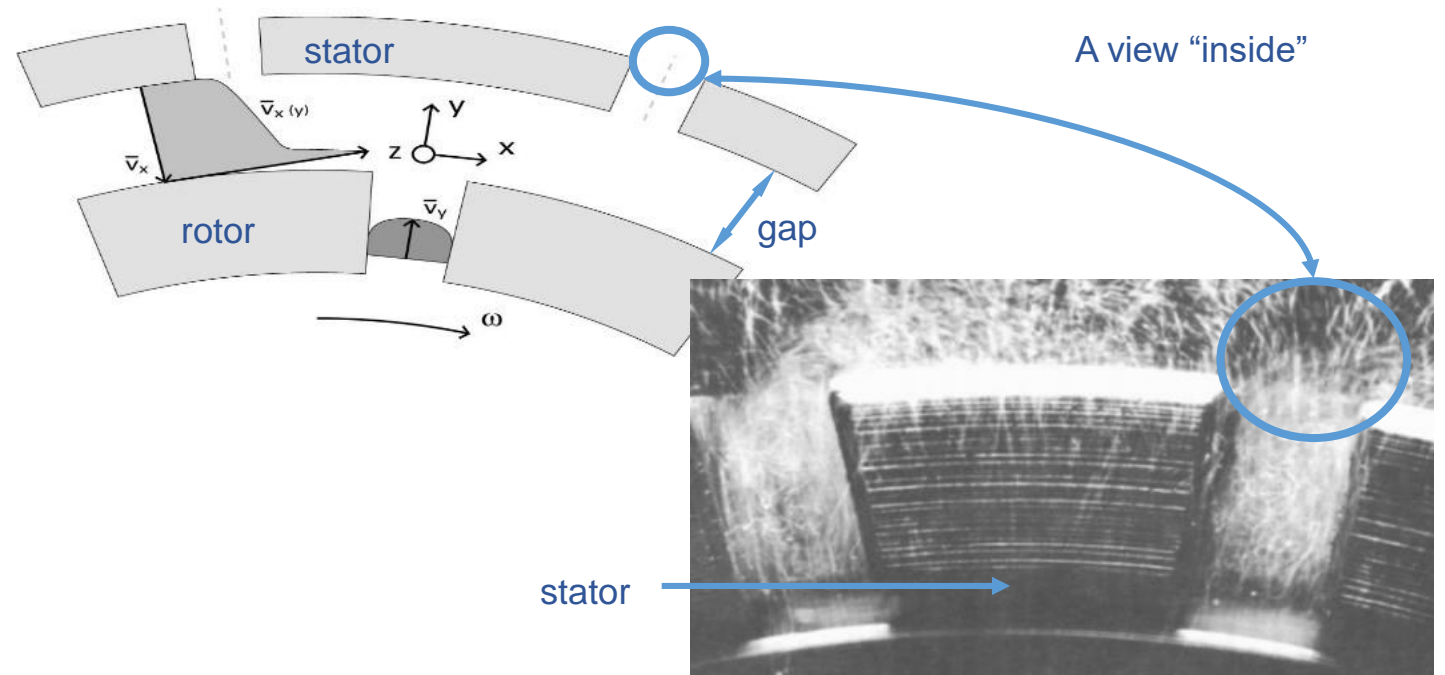


OV5 – ROTOR-STATOR PRINCIPLE

The dispersing head or generator consists of a rotor which rotates at high speed within a stationary stator.

The rotor and the stator have one or more rows of sharp-edged “teeth”.

Because of the high rotor speed, the medium to be processed is **automatically sucked axially into the dispersing head** and is then **pressed radially through the slots of the rotor-stator arrangement**.

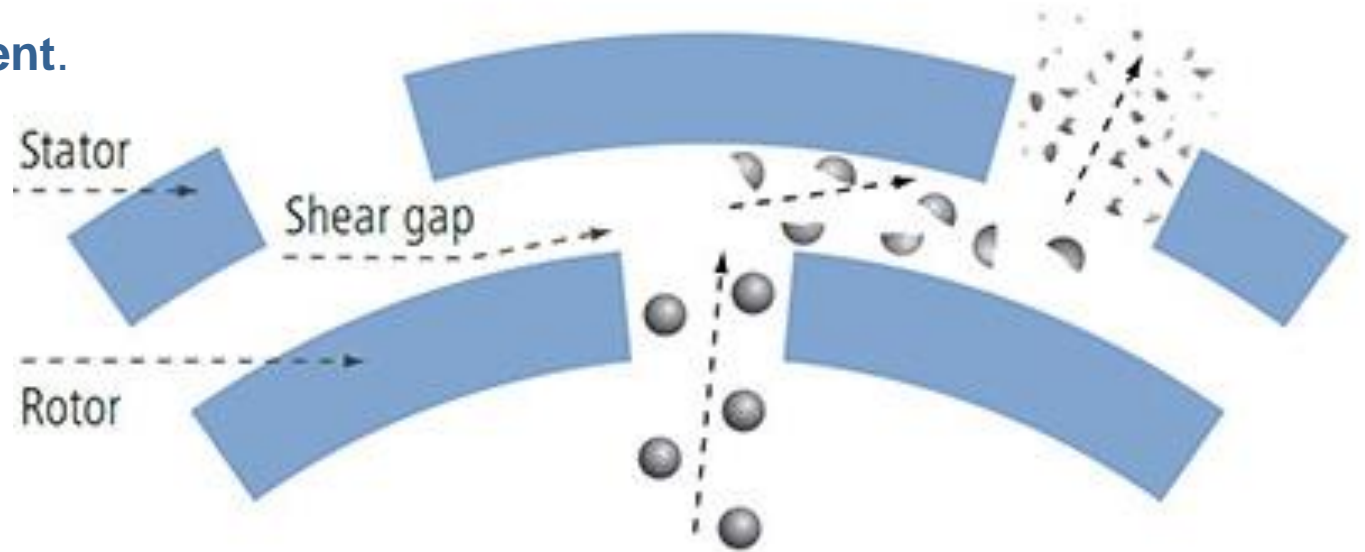


OV5 – ROTOR-STATOR PRINCIPLE

Due to the great forces of acceleration, the material is subjected to **very high shearing** and **thrusting forces**, causing the reduction of the particle size.

The additional high turbulence occurring in the **shearing gap** between the rotor and the stator leads to an **intensive mixing of the material**.

Dispersers are also referred to as “high shear” machines.

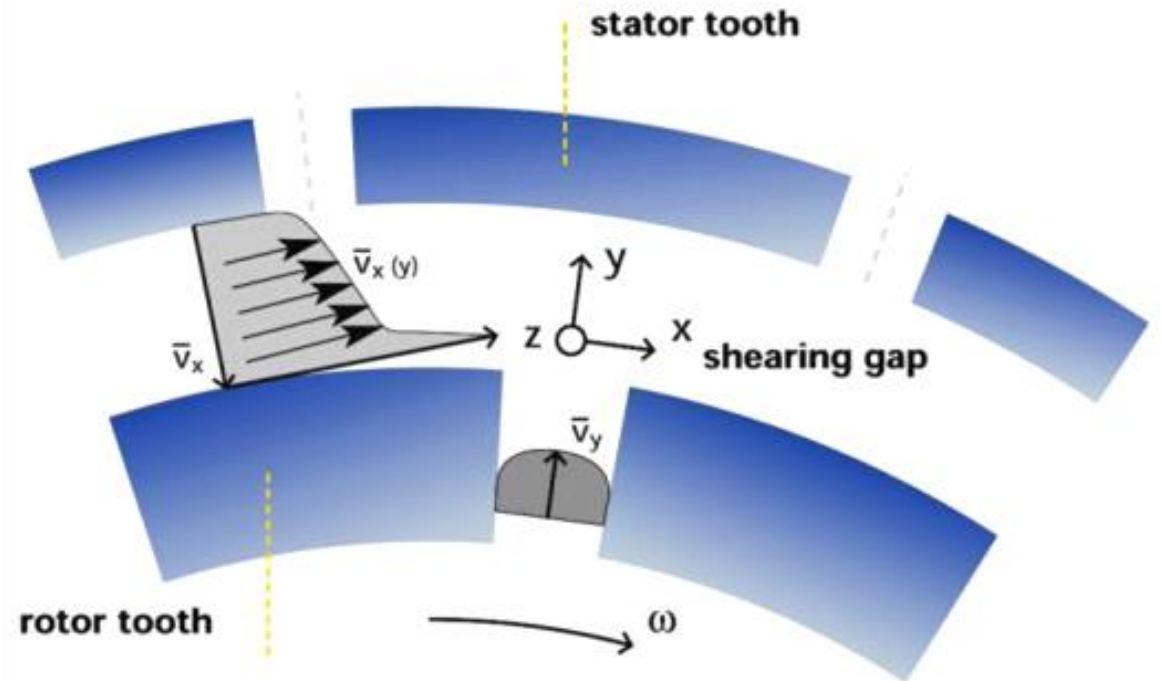

OV5 - SHEARING FORCE


The Shearing force is calculated using the following formula:

$$\text{Shear force} = \frac{\text{Peripheral speed}}{\text{Clearance}}$$

- Peripheral speed is the **speed of the rotor**
- Clearance is the **distance between rotor and stator**. (Usually between 0,3 mm and 1 mm)

For best results, the peripheral speed should be maximised and the clearance minimised.


Shear in Inversely proportional to Viscosity

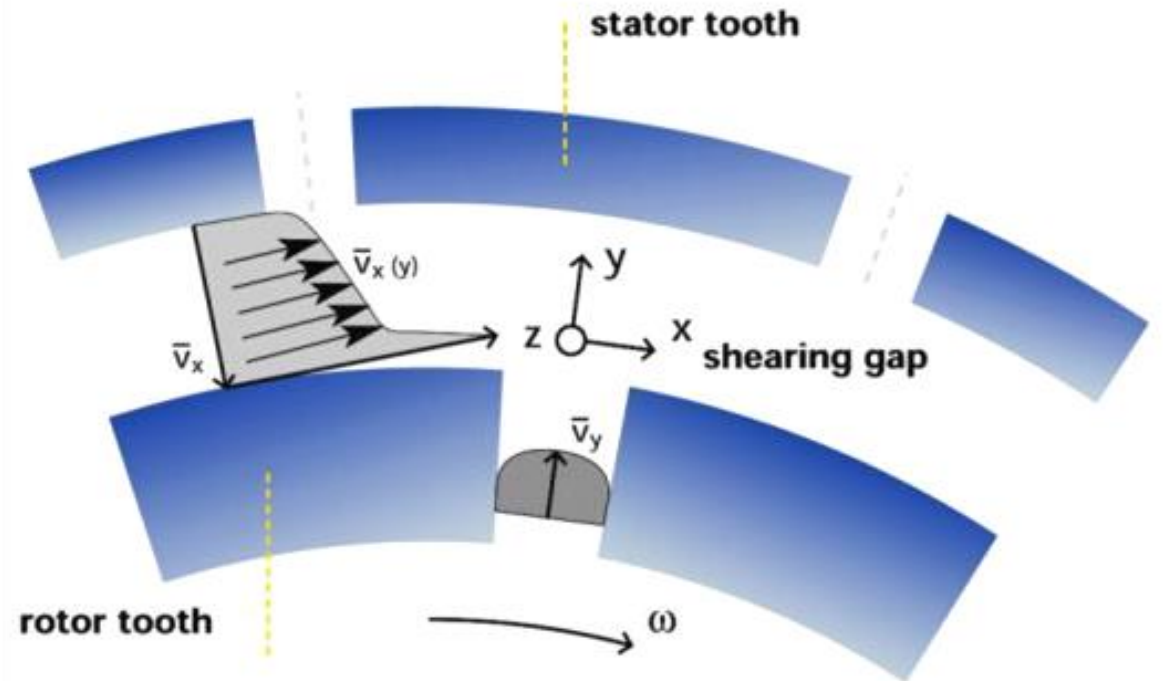

OV5 – ROTOR-STATOR PRINCIPLE

HOW IS FLUID SHEAR INDUCED IN ROTOR/ STATOR?

- A **small portion/parcel** is allowed to enter into the gap (rotor-stator space) before original stream is cut off by the turning rotor. The fluid stream is reduced or braked as the fluid stream “sees” an alternative opening or closing on rotor.
- In the gap, fluid is subjected to **intense 3-dimensional shear** especially at the edge of the rotor teeth. Acceleration and braking - due to presence of teeth and slots.
- Combination of influence from gradient impulse results in shear frequency, shear rate and turbulence

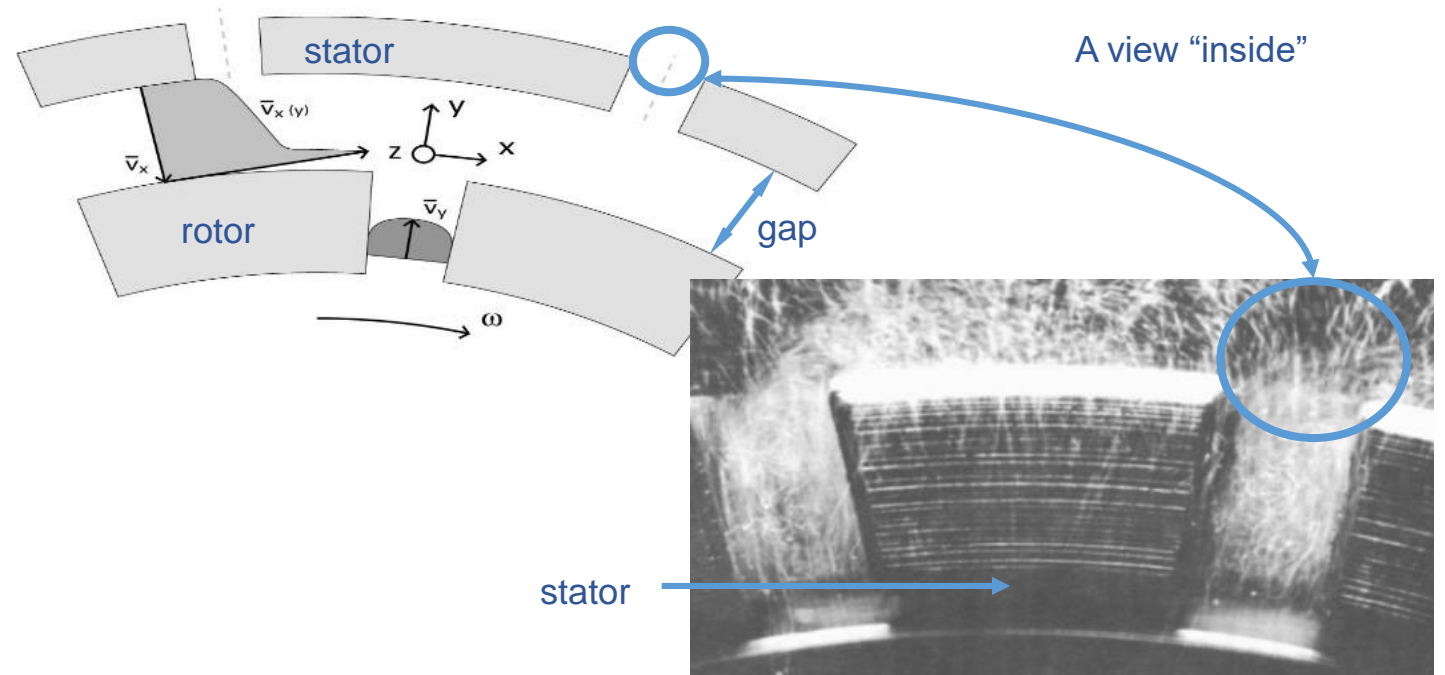

SHEARING FORCE

PERIPHERAL SPEED

Stirrer	0 - 10 m/s
Disperser	10 - 24 m/s

CLEARANCE

Stirrer	0.250 m
(distance from propeller to the vessel wall)	
Disperser	0.005 m

SHEAR FORCE

Stirrer	0-40
Disperser	2000 - 4800



SHEARING FORCE

DISPERSER vs STIRRER

Therefore in the shearing gap of the disperser, its produced about **1000 times the energy is introduced into the medium than when processing by stirring.**

Dispersing times of a few minutes, or even a few seconds, are usually sufficient.

The short processing times is a key factor for not increasing the temperature of the medium due to the energy input (eg. RNA etc..)



OV5 – ROTOR-STATOR PRINCIPLE

BENEFITS

- Time savings
- Energy savings
- Reduction in requirements for additives, stabilizers, emulsifiers

One of the most important parameters for selecting the appropriate instrument is the **volume** to be processed. When selecting an instrument the focus has to be on volume and application

OV5 VOLUMES

- Up to 8 liters as homogenizer
- up to 40 liters as high speed mixer



OV5 – ROTOR-STATOR PRINCIPLE

QUESTION TO ASK

- What is the initial size of the material to be dispersed?
- What is the ultimate fineness required ?
- If the material which is to be dispersed is initially coarse (8-20 mm), it is useful to work with two dispersing tools with different ultimate fineness.

Important !!

Dry running destroys the seal



APPLICATIONS

BIOCHEMISTRY

- Disruption of animal and human tissues, e.g. from eukaryotic cell lines (human/rat/mouse) for 2D gel electrophoresis of cell lysates, enzyme studies, lipid and protein extraction.
- Fragmentation and disruption of muscle and other (vertebrate) tissues (human brain, lymph node, lung and skin); analysis of brain tissue.
- Sample preparation for RNA and DNA extraction; preparation of cell lysates and plasma proteins in vivo

APPLICATIONS

MEDICAL

- Sample preparation of cell dispersion for measurement of all intracellular glutathione from cultured kidney tubule cells;
- Examination of carcinomas from biopsies in cancer research;
- Homogenizing heart, liver, and skeletal tissues down to the subcellular level in the field of experimental surgery.



APPLICATIONS

PHARMACEUTICAL

- Organ extracts
- Tablet / pill coatings
- Serum production
- Herbal extracts
- Gels
- Grinding pills for determination of
- Ointment production
- Color design of end product.



APPLICATIONS

COSMETICS INDUSTRY

- Making emulsions from oil and water
- Productions of lotions
- Hair care articles
- Bath salts
- Sun tan lotions
- Deodorants, gels, nail varnishes lip sticks
- Natural cosmetics,
- Mascara
- Shaving cosmetics
- Tooth paste
- Showering additives, hair dyes, plant extracts, shampoo, skin oils.



APPLICATIONS

FOOD AND FEED INDUSTRY

- Enzymatic treatment
- Preparation of concentrates
- Mixing of beverages
- Stabilizers
- Charging with sugar (e.g. glucose and fructose)
- Milk processing
- Baby food
- Ketchup
- Dietary food
- Feed
- Blending of coloring matter and flavors



APPLICATIONS

CHEMICAL INDUSTRY - PETROCHEMISTRY

- Adding lubrication / lube oil additives
- Sample preparation for determining the water content of crude oil(s)
- Analyzing aging of lubricating oils according to DIN (resistance to oxidation)
- Laboratory-scale mineral and lubricating oil production,
- Production of coated paper, e.g. in the food industry / sector,
- Bitumen for road construction, synthetic bonding agents for making colored asphalts, color coatings or colored isolations.

PAINT AND LACQUER INDUSTRY

- Production of printing links
- Latex and emulsion paint
- Impregnating agents for wood, loading (with) pigments, disintegrating thixotropic materials /matter
- Finest grade homogenizing of water lacquers
- Metal lacquers



OV5 - BENEFITS

- All dispersing tools can be used with the same drive
- Any **rotor/stator can be exchanged** thus no requirement to buy one complete shaft for each application. You buy what really is needed
- **Easy dismantling** and **assembly of shafts** for sterilizing purpose without special tools within a few seconds
- High rpm rates of **10000-30000 rpm**, very high tip speeds result in shorter processing time and finer particle sizes
- Safety switch protects the operator
- Jerk free start and constant rpm rates
- Excellent price/performance ratio

VELP high quality Analytical Dept. at your disposal for fine tuning your application



OV5 - BENEFITS

GENERAL FEATURES

Construction material:	Technopolymer
Very fast assembly and disassembly system for dispersing tools	
Power:	500 W
Weight:	1.3 Kg (2.9 lb)
Dimensions (WxHxD):	70x255x70 mm (2.7x10x2.7 in)

PERFORMANCE

Stirring speed:	from 10'000 to 30'000 rpm
Stirring speed:	graduated scale
Stirring volume max (H2O):	up to 8 liters as homogenizer
Stirring volume max (H2O):	up to 40 liters as high speed mixer
Max viscosity (mPa*s)*:	10'000



OV5 ACCESSORIES



A00000045

H-stand with strap clamp,
bosshead clamp

With this single code is
possible to get all the
necessary tools for
supporting the OV5



A00001300

Support rod and base



A00001300

Double Clamp



A00000044

Strap clamp



A00000279

Stator rotor release kit

**This accessories combined together allow a second configuration for
OV5 support.**

OV5

Example model composition
VSS2CSR2

VS + S2C + SR2



Shaft
VS



Stator
VS2C



Rotor
VSR2



OV5



OV5 – DISPERSING TOOLS



A00000026 – VSS2CSR2
Dispersing Tool for solid/liquid

Applications:
Ceramic Industry – Pharmaceutical Industry –
Petrochemical Industry Sewage Pollution Control



A00000027 – VSS2CCR2
Dispersing Tool with blades for fibrous/stringy materials

Applications:
Ceramic Industry – Pharmaceutical Industry –
Petrochemical Industry Sewage Pollution Control



A00000028 – VSS2CMR2
Dispersing Tool for solid/liquid

Applications:
Ceramic Industry – Paint Industry



A00000029 – VSS2FER2
Dispersing Tool for water, oil, or oil/water emulsions

Applications:
Paper and Tissue – Pharmaceutical Industry – Sewage
pollution control- Paint Industry



A00000031 – VSS2FCR2
Dispersing Tool with blades for fibrous/stringy materials

Applications:
Biotechnology – Paper & Tissue – Food Industry –
Tobacco Industry – Medical – Sewage pollution Control



A00000032 – VSS2FER2
Dispersing Tool solid/liquid media

Applications:
Ceramic Industry – Paper & Tissue – Food Industry –
Cosmetic – Petrochemical – Paint Industry

OV5 – DISPERSING TOOLS



A00000033 – VSS3CSR3
Dispersing Tool solid/liquid media

Applications:
Paper & Tissue – Food Industry – Pharmaceutical –
Medical – Sewage pollution Control



A00000034 – VSS3CCR3
Dispersing Tool with blades for fibrous/stringy materials

Applications:
Paper & Tissue – Food Industry – Pharmaceutical –
Medical – Sewage pollution Control



A00000035 – VSS3CMR3
Dispersing Tool solid/liquid media

Applications:
Ceramic Industry – Paint industry



A00000036 – VSS3CMR2
Stirring shaft

Applications:
Ceramic – Food Industry – Sewage pollution Control



A00000037 – VSS3FER3
Dispersing Tool for water/oil or water/oil emulsions

Applications:
Paper & Tissue – Pharma – Sewage Pollution Control –
Paint Industry



A00000038 – VSS3FSR3
Dispersing Tool for solid/liquid media

Applications:
Paper & Tissue – Pharma – Sewage Pollution Control –
Paint Industry

OV5 – DISPERSING TOOLS



A00000040 – VSS3FMR3

Dispersing Tool for solid/liquid media

Applications:

Ceramic Industry – Food Industry – Cosmetic – Pharma
– Tobacco



A00000041 – VSS4CMR3

Stirring shaft

Applications:

Ceramic Industry – Food Industry – Sewage pollution
control



A00000046 – VSS5CSR4

Dispersing Tool for solid/liquid media

Applications:

Biotechnology – Medical

OV5 – DISPERSING TOOLS

AVAILABLE CONFIGURATIONS

										ULTIMATE FINENESS μm		
i	MODEL	CODE NO	APPLICATION FIELD***	FUNCTION +++	TREATABLE VOLUME (WATER) ml	MAX. CIRCUM. SPEED m/s	Ø ROTOR mm	Ø STATOR mm	TOOL LENGTH mm	MIN/MAX IMMERSION DEPTH mm	SUSPENSION	EMULSION
	*VSS2CSR2	A00000026	CE,IF,PC,SI	A	100-5000	22,7	15	20	220	40/175	10 - 50	1- 10
	VSS2CCR2	A00000027	CT,IA,IT,M,SI	B	100-5000	22,7	15	20	220	40/175	10 - 50	1- 10
	VSS2CMR2	A00000028	CE,VE	A	100-5000	22,7	15	20	220	40/175	10 - 50	1- 10
	*VSS2FER2	A00000029	CT,IF,SI,VE	C	100-5000	22,7	15	20	220	40/175	10 - 50	1- 10
	VSS2FCR2	A00000031	BT,CT,IA,IT,M,SI	B	100-5000	22,7	15	20	220	40/175	10 - 50	1- 10
	VSS2FMR2	A00000032	CE,CT,IA,IC,PC,VE	A	100-5000	22,7	15	20	220	40/175	10 - 50	1- 10
	*VSS3CSR3	A00000033	CT,IA,IF,M,SI	A	100-8000	34,9	23	30	220	40/175	5 - 25	1- 5
	VSS3CCR3	A00000034	CT,IA,IF,M,SI	B	100-8000	34,9	23	30	220	40/175	5 - 25	1- 5
	VSS3CMR3	A00000035	CE,VE	A	100-8000	34,9	23	30	220	40/175	5 - 25	1- 5
	VSS3CMR2	A00000036	CE,IA,SI	D	250-20000	34,9	15	30	220	40/175	High speed mixer	
	*VSS3FER3	A00000037	CT,IF,SI,VE	C	100-8000	34,9	23	30	220	40/175	5 - 25	1- 5
	VSS3FSR3	A00000038	CT,IF,SI,VE	A	100-8000	34,9	23	30	220	40/175	5 - 25	1- 5
	VSS3FMR3	A00000040	CE,IA,IC,IF,IT	A	100-8000	34,9	23	30	220	40/175	5 - 25	1- 5
	*VSS4CMR3	A00000041	CE,IA,SI	D	1000-40000	34,9	23	40	220	40/175	High speed mixer	
	**VSS5CSR4	A00000046	BT,M	A	0,2-50	6,3	4	5	128	10/60	10 - 50	1- 10

* Most used model

** The dispersing tool works with Ø 4 mm rotor and Ø 5 mm stator for microbiological applications (e.g. suitable for Eppendorf, cuvettes,etc.)

*** BT = biotechnology, CE = ceramic industry, CT = paper & tissue industry, IA = food industry, IC = cosmetic industry, IF = pharmaceutical industry, IT = tobacco industry, M = medicine, PC = petrochemistry industry, SI = sewage pollution control, VE = paint industry

**** A = dispersing tool for solid/liquid media, B = dispersing tool with blades for fibrous/stringy materials, C = dispersing tool for water/oil or oil/water emulsions, D = stirring shaft